

CENTRIFUGE **CENTRIC CF48/-R**

TRUE SIMPLICITY - CLASSIC CENTRIFUGES



NEW

MAIN ADVANTAGES:

- > Compact design
- > Operational immediately
- > Ease of use
- > End of run indication

Centrifuges CENTRIC CF 48/-R series, low speed, 4x400 ml, offer an excellent capacity/space ratio and enhanced ergonomics: motorized lid opening, 10 program memories for direct access, AWELight illuminated handle to indicate end of run...

TECHNICAL SPECIFICATIONS		
Model	CENTRIC CF48	CENTRIC CF48-R
Power supply	230V / 50Hz	
Max. speed / RCF	4.500 rpm / 3.820 x g	
Swing-out rotor	4.400 rpm / 3.030 x g	
Angle rotor		
Max. capacity	4 x 400 ml	
Swing-out rotor	32 x 15 ml	
Angle rotor		
Control system	microprocessor	
Speed range	500 to 4.500 rpm	
Speed precision	20 rpm	
Temperature control range	-	-20°C to +60°C
Time	0 - 99h50 and continuous	
Acceleration / Deceleration	10 / 10	
Program settings	10 programs	
Lid opening	motorized	
Motor type	induction	
Certification	EC marked: IEC 1010-1; IEC 1010-2-020 directives	
Dimensions (HxWxD)	280 x 420 x 520 mm	280 x 420 x 750 mm
Crating (HxWxD)	460 x 610 x 650 mm	460 x 610 x 890 mm
Weight (Net/Crated)	50 / 56 kg	70 / 79 kg



 **SMART LIGHT**
End of run light indication

 **SMART FREEZE**
Superior cooling cycles

 **SMART OPEN**
Motorized lid opening

 **SMART PURE**
High stability, low vibrations for better results

 **SMART LOOK**
Compact design

 **SMART EASY**
Intuitive programming

 **SMART LOCK**
Rotor exchange without tools


 **SMART SEAL**
One-handed biosafety lid closure

Distribuito da:
Zetalab s.r.l.

Via Umberto Giordano, 5 - 35132 Padova
Tel 049 2021144 - Fax 049 2021143
www.zetalab.it - email: info@zetalab.it

Zetalab.it

A WIDE RANGE TO SUIT ALL YOUR APPLICATIONS

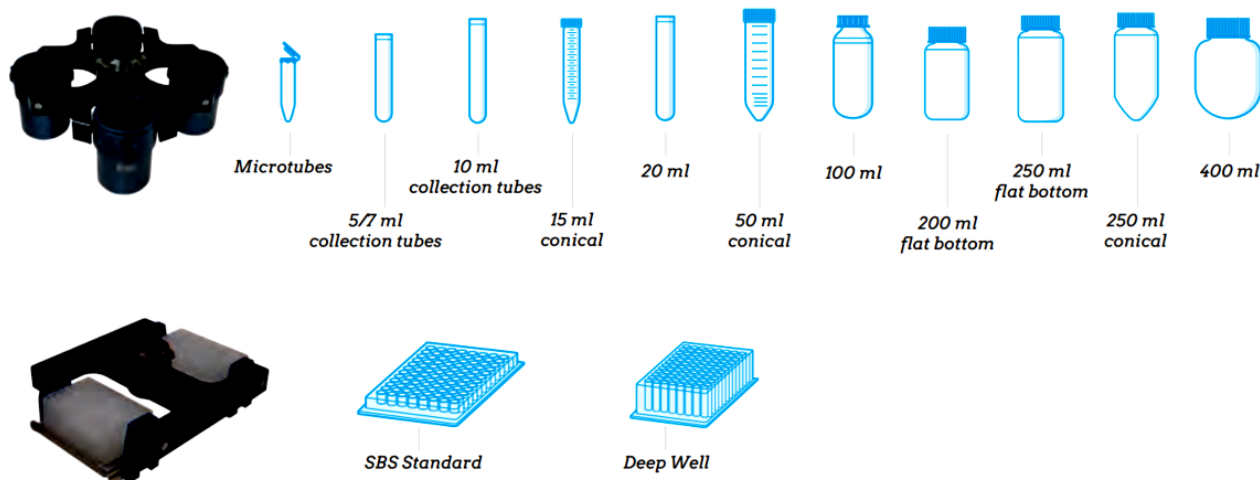

Swing-out rotors		
SL-400	4 place swing-out rotor	4 x 400 ml
N-400 S	Buckets for SL-400	Set of 4
BL-400 S	Lid for N-400 S	Set of 4
	Inserts for 13 x 1.5/2 ml microtubes (52 microtubes / rotor)	Set of 4
	Inserts for 14 x 5-7 ml collection tubes (56 collection tubes / rotor)	Set of 4
	Inserts for 15 x 5-7 ml BD collection tubes (60 collection tubes / rotor)	Set of 4
	Inserts for 19 x 5-7 ml Terumo coll. tubes (76 Terumo collection tubes / rotor)	Set of 4
	Inserts for 12 x 10 ml collection tubes (48 collection tubes / rotor)	Set of 4
	Inserts for 7 x 15 ml conical tubes (28 conical tubes / rotor)	Set of 4
	Inserts for 5 x 20 ml tubes (20 tubes / rotor)	Set of 4
	Inserts for 4 x 50 ml conical tube (16 conical tubes / rotor)	Set of 4
	Inserts for 1 x 100 ml tube dia. 52 mm (4 tubes dia. 52 mm / rotor)	Set of 4
	Inserts for 1 x 250 ml tube (4 tubes / rotor)	Set of 4
SL 48-MTP	Swing-out rotor for microtitration buckets	
	2 Microtitration buckets (2 x 3 microplates / rotor)	Set of 2